



How to use **Awesome Floss**[®] dental floss*

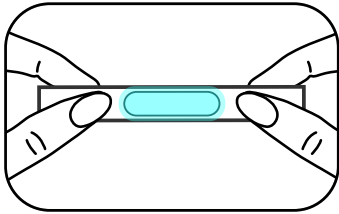


Step 1 Squeeze one end of **Awesome Floss**[®] dental strip firmly between the thumb and index finger of one hand. Squeeze the other end of the **Awesome Floss**[®] dental strip firmly between your other thumb and index finger so that your thumbs are about 1¼ or 1½ inches (3 – 4 cm) apart.


Note: Do *NOT* directly touch or pre-wet the colorful see-thru film *in the middle* of the strip and on *both sides* of the strip. The *breath-freshening colorful center will dissolve* when guided between your teeth revealing *two strands* of floss, one strand for the upper teeth and one for the lower teeth.

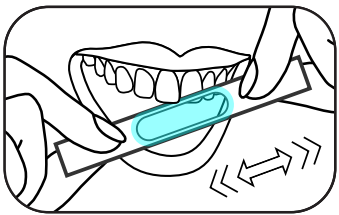


Step 2 Glide the middle of just one strand back and forth between two adjacent teeth.

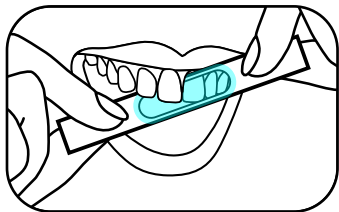


Step 3 Ease the *one strand* between teeth *until one edge* of that strand *gently touches below your gum-line*.

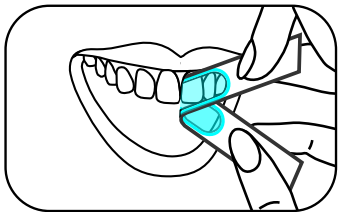

Step 4 Hug that strand tightly against one tooth.

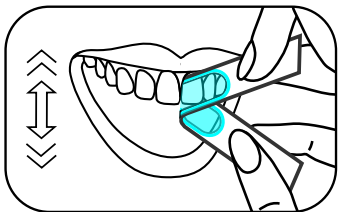

Step 5 Gently move the *hugging strand up and down* (and optionally back and forth) several times to *properly clean above and below your gum-line*.

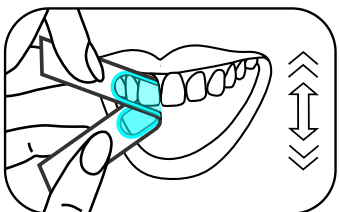

Step 6 Repeat Steps 4 and 5 on the adjacent tooth.

Step 7 Repeat Steps 2 thru 6 for every pair of contacting or adjacent teeth. To clean *behind* each back tooth simply repeat Steps 4 and 5.